
TissUUmaps

Release 3.0

Nicolas Pielawski	Axel Andersson
Christophe Avenel	Andrea Behanova
Eduard Chelebian	Anna Klemm
Leslie Solorzano	Fredrik Nysjö
	Carolina Wählby

Apr 14, 2022

CONTENTS:

1	Introduction	2
1.1	About TissUMaps	2
1.2	Installation	2
1.3	Citing TissUMaps	3
1.4	Changelog	3
2	Getting started	6
2.1	Images	6
2.2	Markers	6
2.3	Regions	7
2.4	Projects	7
2.5	Exporting screenshots	14
2.6	Plugins	15
3	Sharing projects	17
3.1	Apache server	17
3.2	Docker container	17
4	Advanced usage	18
4.1	Jupyter notebooks	18
4.2	Napari	19
4.3	AnnData	20
	Index	21

This page hosts the documentation for TissUMaps 3.0. You can find a pdf version of this documentation [here](#).

For more information on the TissUMaps project, including video tutorials and demos, visit our website: <https://tissuums.github.io>.

Work in progress!

This page is mostly empty for now. We are working actively on writing this documentation, more content will be available soon!

INTRODUCTION

1.1 About TissUMaps

[TissUMaps](#) is a free and open source browser-based tool for GPU-accelerated visualization and interactive exploration of tens of millions of datapoints overlaying tissue samples. Users can visualize markers and regions, explore spatial statistics and quantitative analyses of tissue morphology, and assess the quality of decoding in situ transcriptomics data. TissUMaps provides instant multi-resolution image viewing, can be customized, shared, and also integrated in Jupyter Notebooks. We envision TissUMaps to contribute to broader dissemination and flexible sharing of large-scale spatial omics data.

Currently, microscopy data can be cumbersome to share: physically transferring the images is often necessary and dedicated software must be installed. Instead, researchers can now share their findings with a simple link to a website running TissUMaps. The images are loaded in real time, together with annotations, markers, and masks that may also be modified by the user. We also provide tools for quality control and image processing. The software is designed to display and interact with images at multiple resolutions and large numbers of markers, especially data from spatially resolved omics techniques and tissue atlases. TissUMaps is compatible with many different bioimage informatics tools, and provides new ways to develop insights when exploring and sharing data.

You can access the [TissUMaps project gallery](#) with interactive examples to explore data from in situ sequencing and spatial transcriptomics experiments and view localized quantification of cell and tissue morphology, including links to publications. For seeing examples of TissUMaps compatibility with other platforms you can access the [tutorials page](#).

1.2 Installation

[TissUMaps](#) is a browser-based tool for fast visualization and exploration of millions of data points overlaying a tissue sample. TissUMaps can be used as a web service or locally in your computer, and allows users to share regions of interest and local statistics.

1.2.1 Windows installation

1. Download the Windows Installer from [the last release](#) and install it. Note that the installer is not signed yet and may trigger warnings from the browser and from the firewall. You can safely pass these warnings.

1.2.2 PIP installation (for Linux and Mac)

1. Install libvips for your system: <https://www.libvips.org/install.html>

An easy way to install libvips is to use an [Anaconda](#) environment with libvips:

```
conda create -y -n tissuumaps_env -c conda-forge python=3.9 libvips
conda activate tissuumaps_env
```

2. Install the TissUMaps library using pip:

```
pip install "TissUMaps[full]"
```

3. Start the TissUMaps user interface:

```
tissuumaps
```

4. Or start TissUMaps as a local server:

```
tissuumaps_server path_to_your_images
```

And open <http://127.0.0.1:5000/> in your favorite browser.

1.3 Citing TissUMaps

Please cite our [preprint](#) on bioRxiv if using TissUMaps in your work:

TissUMaps 3: Interactive visualization and quality assessment of large-scale spatial omics data. *Nicolas Pielawski, Axel Andersson, Christophe Avenel, Andrea Behanova, Eduard Chelebian, Anna Klemm, Fredrik Nysjö, Leslie Solorzano, Carolina Wählby*, bioRxiv 2022.01.28.478131; doi: <https://doi.org/10.1101/2022.01.28.478131>.

1.4 Changelog

1.4.1 3.0.8.5

- Minor fixes.

1.4.2 3.0.8.4

- Add tiling to viewport capture for higher resolution output
- Increase resolution of markers on high resolution devices
- Fix jumps on pan with mouse gesture (mobile)
- Add fix for bright image canvas on Safari
- Add an option to remove markers' outlines.

1.4.3 3.0.8.3

- Fix png artifact in Firefox, by generating jpg tiles.

1.4.4 3.0.8.2

- Add high resolution capture of viewport, up to 4096x4096 pixels.

1.4.5 3.0.8.1

- Fix multiple dataset alignment when no background image

1.4.6 3.0.8

- Fix black images generated by VIPS
- Fix Linux and Mac open of captures
- Auto save datasets as buttons when saving tmap projects
- Add mpp (microns per pixel) option in tmap files, to add scale bar to viewer
- Make region line thickness depend on zoom level
- Add compatibility with JupyterLab
- Add opacity per marker option

1.4.7 3.0.7

- Add menu to load plugins through an update-site

1.4.8 3.0.6

- Fix multiple plugins opening always last plugin
- Move to OpenSeadragon 3.0.0
- Add tooltip format in Advanced Settings
- Add drag and drop to open CSV files and images
- Add “Add layer” button for flask version
- Add viewport capture

1.4.9 3.0.5

- Move csv loading to Papa Parse streaming, to allow better memory management

1.4.10 3.0.4

- Add filtering of markers

1.4.11 3.0

- Add tissuumaps.jupyter module

GETTING STARTED

2.1 Images

2.1.1 Supported image formats

TissUUmapi can read whole slide images in any format recognized by the OpenSlide library:

- Aperio (.svs, .tif)
- Hamamatsu (.ndpi, .vms, .vmu)
- Leica (.scn)
- MIRAX (.mrxs)
- Philips (.tiff)
- Sakura (.svslide)
- Trestle (.tif)
- Ventana (.bif, .tif)
- Generic tiled TIFF (.tif)

TissUUmapi will automatically convert any other format into a pyramidal tiff (in a temporary .tissuumapi folder created in the original image folder) using vips.

If your image fails to open, try converting it to `tif` format using an external tool.

2.1.2 Load images

2.1.3 Apply filters

2.2 Markers

2.2.1 Supported marker format

TissUUmapi can read CSV (Comma Separated Values) files with a header row, and at least spatial coordinate columns (X and Y). CSV files are not limited in the number of columns or number of rows. Other columns can contain information for displaying markers (key to group markers, color, size, shape, piecharts, etc.)

CSV files can be exported from any spreadsheet program, or any programming language (Python, R, etc.)

2.2.2 Load markers

2.2.3 Markers settings

File and coordinates

Render options

Advanced options

Table of markers

2.3 Regions

2.3.1 Supported region formats

TissUMaps can read and write region files in the [GeoJSON](#) format.

Only a subset of the GeoJSON format is supported, as TissUMaps uses only polygonal regions:

Main types:

- Feature
- FeatureCollection
- GeometryCollection

Geometries:

- Polygon
- Multipolygon

The coordinate system must be the same as the image and marker coordinate systems.

2.3.2 Draw Regions

2.3.3 Analyze Regions

2.3.4 Load Regions

2.3.5 Export Regions

2.4 Projects

2.4.1 Saving and loading projects

2.4.2 The tmap file format

The tmap format contains image layers, saved markers, regions, and settings. It is highly recommended to create tmap files by saving projects from TissUMaps applications, but you can also edit the files manually to add or change project's settings.

The tmap format uses json, with the following specifications:

Tmap specifications

type	<i>object</i>		
properties			
• filename	Name of the project		
	type	<i>string</i>	
• layers	type	<i>array</i>	
	items		
	•	<i>Layer</i>	
• layerOpacities	type	<i>object</i>	
	patternProperties		
	• ^[0-9]+\$	type	<i>integer</i>
• layerVisibilities	type	<i>object</i>	
	patternProperties		
	• ^[0-9]+\$	type	<i>boolean</i>
• layerFilters	type	<i>object</i>	
	patternProperties		
	• ^[0-9]+\$	<i>LayerFilter</i>	
• filters	type	<i>array</i>	
	items		
	•	<i>Filter</i>	
• compositeMode	type	<i>string</i>	
• markerFiles	type	<i>array</i>	
	items		
	•	<i>MarkerFile</i>	
• regions	GeoJSON object, see <i>Regions section</i> .		
	type	<i>object</i>	
• regionFile	type	<i>string</i>	
• regionFiles	type	<i>array</i>	
	items		
• plugins	type	<i>array</i>	
	items		

continues on next page

Table 1 – continued from previous page

	•	type	<i>string</i>
• hideTabs	Hide tabs of markers dataset. Only use when you have a unique marker tab.		
	type	<i>boolean</i>	
• settings	type	<i>array</i>	
	items		
	•	<i>Setting</i>	

Layer

type	<i>object</i>		
properties			
• name	type	<i>string</i>	
• tileSource	type	<i>string</i>	

LayerFilter

type	<i>array</i>		
items			
•	type	<i>object</i>	
	properties		
	• name	<i>Filter</i>	
	• value	type	<i>string</i>

Filter

enum	Color, Brightness, Exposure, Hue, Contrast, Vibrance, Noise, Saturation, Gamma, Invert, Greyscale, Threshold, Erosion, Dilation
------	---

MarkerFile

Description of settings and GUI objects for a marker dataset loaded from CSV file		
type	object	
properties		
• title	Name of marker button	
	type	string
• comment	Optional description text shown next to marker button	
	type	string
	default	
• name	Name of marker tab	
	type	string
• autoLoad	If the CSV file for the marker dataset should be automatically loaded when the TMAP project is opened. If this is false, the user has to instead click on the marker button to load the dataset.	
	type	boolean
	default	false
• uid	A unique identifier used internally for the marker dataset	
	type	string
• expectedHeader	ExpectedHeader	
• expectedRadios	ExpectedRadios	
• path	Relative file path to CSV file for marker dataset	
	type	string
• settings	type	array
	items	
	•	Setting

ExpectedHeader

Input field values for settings in a marker tab		
type	object	
properties		
• X	Name of CSV column used for X-coordinate	
	type	string
• Y	Name of CSV column used for Y-coordinate	
	type	string
• gb_col	Name of CSV column to group markers by	
	type	string

continues on next page

Table 2 – continued from previous page

	default	null
• gb_name	Name of alternative CSV column to show for group instead of gb_col	
	type	<i>string</i>
	default	null
• cb_cmap	Name of colorscale to use for color mapping	
	type	<i>string</i>
	default	
• cb_col	Name of CSV column containing scalar values for color mapping or hexadecimal colors	
	type	<i>string</i>
	default	null
• cb_gr_dict	TODO	
	type	<i>string</i>
	default	
• scale_col	Name of CSV column containing scalar values for changing the size of markers	
	type	<i>string</i>
	default	null
• scale_factor	Numerical value for a fixed scale factor to be applied to markers	
	type	<i>string</i>
	default	1
• pie_col	Name of CSV column containing data for pie chart sectors	
	type	<i>string</i>
	default	null
• pie_dict	TODO	
	type	<i>string</i>
	default	
• shape_col	Name of CSV column containing a name or an index for marker shape	
	type	<i>string</i>
	default	null
• shape_fixed	Name of a fixed marker shape to be use for all markers	
	type	<i>string</i>
	default	cross
• shape_gr_dict	TODO	
	type	<i>string</i>
	default	
• opacity_col	Name of CSV column containing scalar values for opacities	
	type	<i>string</i>
	default	null
• opacity	Numerical value for a fixed opacity factor to be applied to markers	
	type	<i>string</i>
	default	1
• tooltip_fmt	TODO	
	type	<i>string</i>
	default	

ExpectedRadios

Radio button state and checkbox state for settings in a marker tab		
type	<i>object</i>	
properties		
• cb_col	type	<i>boolean</i>
	default	false
• cb_gr	type	<i>boolean</i>
	default	true
• cb_gr_rand	type	<i>boolean</i>
	default	false
• cb_gr_dict	type	<i>boolean</i>
	default	false
• cb_gr_key	type	<i>boolean</i>
	default	true
• pie_check	type	<i>boolean</i>
	default	false
• scale_check	type	<i>boolean</i>
	default	false
• shape_gr	type	<i>boolean</i>
	default	true
• shape_gr_rand	type	<i>boolean</i>
	default	true
• shape_gr_dict	type	<i>boolean</i>
	default	false
• shape_col	type	<i>boolean</i>
	default	false
• shape_fixed	type	<i>boolean</i>
	default	false
• opacity_check	type	<i>boolean</i>
	default	false
• _no_outline	type	<i>boolean</i>
	default	false

Setting

type	<i>object</i>	
properties		
• function	type	<i>string</i>
• module	type	<i>string</i>
• value	type	<i>number</i>

Example of tmap file

```
{
  "filename": "TissUMaps_Example.tmap",
  "layers": [
    {
      "name": "Round1_A.tif",
      "tileSource": "images/Round1_A.tif.dzi"
    },
    {
      "name": "Round1_C.tif",
      "tileSource": "images/Round1_C.tif.dzi"
    }
  ],
  "layerOpacities": {
    "0": "1",
    "1": "1"
  },
  "layerVisibilities": {
    "0": true,
    "1": false,
  },
  "layerFilters": {
    "0": [
      {
        "name": "Color",
        "value": "0,100,0"
      }
    ],
    "1": [
      {
        "name": "Color",
        "value": "0,100,0"
      }
    ]
  },
  "filters": [
    "Color"
  ],
  "compositeMode": "lighter",
  "markerFiles": [
    {
      "autoLoad": false,
      "comment": "",
      "expectedHeader": {
        "X": "global_x",
        "Y": "global_y",
        "cb_cmap": "",
        "cb_col": "null",
        "cb_gr_dict": "",
        "gb_col": "Gene",
        "gb_name": "",
        "opacity": "1",

```

(continues on next page)

(continued from previous page)

```

        "opacity_col": "null",
        "pie_col": "null",
        "pie_dict": "",
        "scale_col": "null",
        "scale_factor": "0.5",
        "shape_col": "null",
        "shape_fixed": "cross",
        "shape_gr_dict": "",
        "tooltip_fmt": ""
    },
    "expectedRadios": {
        "cb_col": false,
        "cb_gr": true,
        "cb_gr_dict": false,
        "cb_gr_key": true,
        "cb_gr_rand": false,
        "pie_check": false,
        "scale_check": false,
        "shape_col": false,
        "shape_fixed": false,
        "shape_gr": true,
        "shape_gr_dict": false,
        "shape_gr_rand": true,
        "opacity_check": false
    },
    "name": " markers",
    "path": "./istdeco_codes_n.csv",
    "title": "Download markers",
    "uid": "uniquetab"
    }
],
"regions": {},
"plugins": [
    "Spot_Inspector"
],
"hideTabs": true,
"settings": []
}

```

2.5 Exporting screenshots

TissUUmapi allows high resolution capture of the image viewport. Go to Menu > File > Capture viewport and chose a zoom factor for export (1 = screen resolution).

The screen capture will contain all filtered layers, markers, and regions. Note that legends will not be part of the export and must be added manually.

2.6 Plugins


2.6.1 Load plugins

2.6.2 Make your own plugin

Download the Plugin Template python and javascript files from the [Plugin Update Site](#) and put both files in your local folder \$USER_PATH/.tissuumaps/plugins/. You can then change the plugin name and add your own options and functions.

Javascript file

When loading a plugin, the function `PluginName.init(container)` will be called. The `container` is an html Element that will be added to the plugin menu. Use this element to add options and texts related to your plugin.

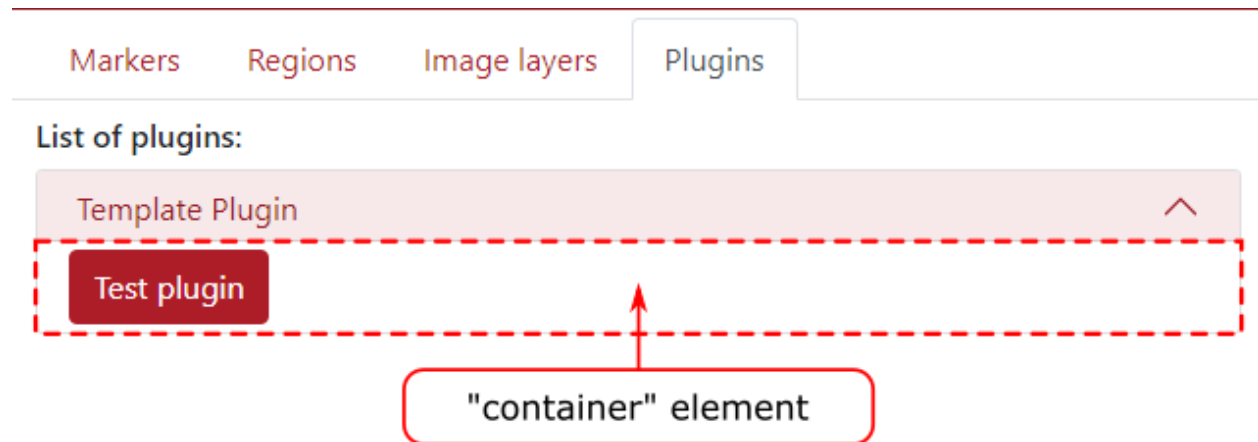

Here is a minimal example of plugin:

```
var Plugin_template;
Plugin_template = {
  name:"Template Plugin"
}

/**
 * This method is called when the document is loaded.
 * The container element is a div where the plugin options will be displayed. */
Plugin_template.init = function (container) {
  container.innerHTML = "Hello world";
}
```

You can access the TissUUmapi javascript API [here](#).

Python file

You only need to use the Python file if your plugin needs to do processing on the server side. For pure javascript plugins, you can leave this file empty.

The python file should implement the class `Plugin`:

```
class Plugin():
    def __init__(self, app):
        self.app = app
```

The `app` object being the flask application running the TissUMaps server.

You can call a Python method inside the `Plugin` class from Javascript using Ajax and the Python API. The endpoint for a method `methodName` of the plugin `PluginName` will be: `/plugins/pluginName/functionName`. Data can be transmitted through Ajax as stringified JSON, and will be available as a parameter inside the method.

See the [Plugin Template](#) for a working example of Javascript / Python communication.

SHARING PROJECTS

3.1 Apache server

TissUMaps projects can be exported into static webpages, that can be uploaded to any Apache server.

1. Save your project from TissUMaps (menu > File > Save project)
2. Export to static page (menu > File > Export to static webpage)
3. Copy the exported folder on your Apache server

3.2 Docker container

1. Start the docker container `cavenel/tissuums:latest` from Docker Hub:

```
docker run -it -p 56733:80 --name=tissuums -v /path/to/local/images:/mnt/data cavenel/  
↪tissuums:latest
```

1. Place your images in the local folder `/path/to/local/images/share`.
2. Open <http://127.0.0.1:56733/> in your favorite browser.

ADVANCED USAGE

4.1 Jupyter notebooks

TissUMaps can easily be used inside a Jupyter Notebook or Jupyter Lab.

Simple example to load an image in TissUMaps:

```
import tissuumaps.jupyter as tj
viewer = tj.loaddata(["image.png"])

viewer.screenshot()
```

4.1.1 tissuumaps.jupyter

Module used to run TissUMaps from a Jupyter Notebook or from Jupyter Lab.

`tissuumaps.jupyter.opentmap(path, port=5100, host='localhost', height=700)`

Open a tmap project

Parameters

- **path** (*str*) – The path to a tmap file
- **port** (*int*) – The port to run the TissUMaps server
- **host** (*str*) – The host to run the TissUMaps server
- **height** (*int*) – The height of the jupyter iframe

Returns The TissUMaps viewer

Return type *TissUMapsViewer*

```
tissuumaps.jupyter.loaddata(images=[], csvFiles=[], xSelector='x', ySelector='y', keySelector=None,
                             nameSelector=None, colorSelector=None, piechartSelector=None,
                             shapeSelector=None, scaleSelector=None, fixedShape=None, scaleFactor=1,
                             colormap=None, compositeMode='source-over', boundingBox=None,
                             port=5100, host='localhost', height=700, tmapFilename='_project',
                             plugins=[])
```

Load data in TissUMaps

Parameters

- **images** (*list* | *str*) – List of images or single image to display
- **csvFiles** (*list* | *str*) – List of csv files or single csv file to display

- **xSelector** (*str*) – Name of the csv column defining the X coordinates
- **ySelector** (*str*) – Name of the csv column defining the Y coordinates
- **keySelector** (*str*) – Name of the csv column defining the grouping key
- **nameSelector** (*str*) – Name of the csv column defining the group name
- **colorSelector** (*str*) – Name of the csv column defining the group color
- **piechartSelector** (*str*) – Name of the csv column defining pie-charts
- **shapeSelector** (*str*) – Name of the csv column defining markers' shape
- **scaleSelector** (*str*) – Name of the csv column defining markers' scale
- **fixedShape** (*int*) – Name of the markers' shape
- **scaleFactor** (*int*) – Global scale of markers
- **colormap** (*str*) – Name of the colormap used if colorSelector is set
- **compositeMode** – (*str*): Composite mode used for images
- **boundingBox** (*list*) – [X,Y,W,H] of the bounding box to display
- **port** (*int*) – The port to run the TissUMaps server
- **host** (*str*) – The host to run the TissUMaps server
- **height** (*int*) – The height of the jupyter iframe
- **tmapFilename** (*str*) – Name of the project file that will be created
- **plugins** (*list*) – List of plugins to add to the tmap project

Returns The TissUMaps viewer

Return type *TissUMapsViewer*

class tissumaps.jupyter.**TissUMapsViewer**(*server, image, height=700*)

Class representing a TissUMaps viewer instance

screenshot()

Capture the TissUMaps viewport and display image in the Notebook.

class tissumaps.jupyter.**TissUMapsServer**(*slideDir, port=5000, host='0.0.0.0'*)

Class representing a TissUMaps server instance

4.2 Napari

Napari features an important hub containing 118 plugins at the time of writing, many of them expanding further the capabilities of Napari when dealing with biomedical imaging. We thus created our own plugin to allow users to work in Napari, benefit from the tools, scripting and existing plugins, and easily visualize and share the output of their research through TissUMaps.

The [Napari-TissUMaps plugin](#) is available on Napari Hub which makes the installation trivial: from the Napari install/uninstall plugins menu, the `napari-tissumaps` appears in the list and can be installed with a single click. Alternatively, the plugin can be installed with the Python package manager: `pip install napari-tissumaps`.

The plugin can export all standard Napari layers, such as images, labels, points, and shapes and preserves the metadata (opacity, visibility), but also the objects parameters (e.g.: label colors, marker colors and symbols, etc...). To export a TissUMaps project, care must be taken to save all layers of interest and type in a name with the extension `.tmap`, e.g.: `myProject.tmap`. This is important for Napari to delegate the saving of the files to the plugin. A folder is created

and contains all the necessary files and can be loaded in the TissUMaps server, software, Jupyter Notebook, or shared with the community.

The project folders generated by the plugin contain the metadata in a `main.tmap` file, along with folders for each Napari layer types: images, labels, points and regions. Images and labels are saved as plain tif images, points are saved as CSV files, and shapes are saved as GeoJSON. We hope that the use of a simple structure and widespread file formats can simplify the modifying and updating of the TissUMaps project when prototyping with e.g. Jupyter Notebooks. The source code is available at <https://github.com/TissUMaps/napari-tissuums> under the permissive MIT license. A demonstration of the Cellpose plugin of Napari being exported to the TissUMaps web viewer is available at: <https://tissuums.github.io/tutorials/#napari>.

4.3 AnnData

Work in progress

INDEX

L

`loaddata()` (*in module `tissuumaps.jupyter`*), 18

M

module

`tissuumaps.jupyter`, 18

O

`opentmap()` (*in module `tissuumaps.jupyter`*), 18

S

`screenshot()` (*`tissuumaps.jupyter.TissUUmapiViewer` method*), 19

T

`tissuumaps.jupyter`

module, 18

`TissUUmapiServer` (*class in `tissuumaps.jupyter`*), 19

`TissUUmapiViewer` (*class in `tissuumaps.jupyter`*), 19